




Electric Examination Table

GANAS 2

We Provide Customized Overall Solutions
for Medical Equipment

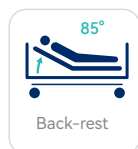

PRODUCT SPECIFICATIONS

External Size (LxW)	1890x680 mm
Height Adjustment	470-910 mm
Head up/down	85°/40°
Leg up/down	48°/45°
Bearing Capacity	150 kg
Headrest Hole	150 mm

Input Voltage	110-230 V
Output Power	24 W

Pillow
Motor
Gas Spring
3" Casters
Brake Pedal
Hand controller
Manual Operator

1pc
1pc
2pcs
4pcs
4pcs
1pc
1pc

Bed Surface Foot Lift
Roll Up Tissue Holder

Optional
Optional



BED SURFACE

Two-fold, ergonomic design; adjustable height for patient convenience, personalized design aids diagnosis.



DUAL USE


Functions as both bed and chair; high-bearing, durable stability; safe for heavy pressure impact.



Accessories

Bed Frame

The bed frame surface is painted with 1.5mm powdercoated steel, and welded by Panasonic Robotic seamlessly. Strong corrosion resistance, environmentally-friendly, bright surface, maintain an overall beautiful appearance. With this structure, the whole bed maintains more stable and makes patients feel more safe.

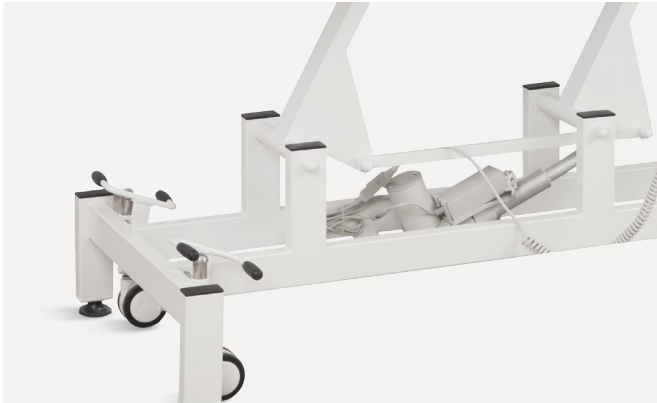

Mattress

The bed platform is integrally made of the high-end composite leather and hand-made cotton mattress, of which the curve design is ergonomic, elastic and soft, breathable, environmentally-friendly, easy to be cleaned, anti-bacterial, impact-resistant.



Pneumatic Pole support

Two air springs are equipped on the back and legs to provide users with a quiet and safe environment during diagnosis and avoid interference.



Control System

CE certified linear actuator motor, one-button free lift, low noise, and long service life. Manual (foot control) electric overall lifting, manual back-rest lifting, leg-rest lifting adjustment, height adjustment.



Hand Controller

The bed is equipped with Hand-held controller, so that the optimal height is easy to be adjusted by pressing the single button gently. The hook on the back is convenient to take down for patients and nurses.



Bed foot equipment

Equipped with 4pcs 3-inch double-sided wheels, 2 function side brake with 4 pedals, which are convenient to move and fix. The design of 360° rotating pulley made of polyurethane PU is strong, durable, non-slip, strong in pressure resistance, and good in stability.



Assorted Pillow

The neoteric design of the pillow is narrow on the top and wide on the bottom. Its well-shaped and elegant design matches with modern aesthetic ideas and can blend into the environment harmoniously. Pillow can be disassembled according to demand.



Round Covered Edge

The bed platform was covered in 360°. The soft-touch of surface would be safer for users.



Design

Bed surface two fold three ergonomic design, adjust the height according to the patient's height, personalized design help diagnosis. Dual use of bed and chair, high bearing lasting stability, no fear of heavy pressure impact, to ensure safe use.



Comprehensive Hospital Design Solutions




Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

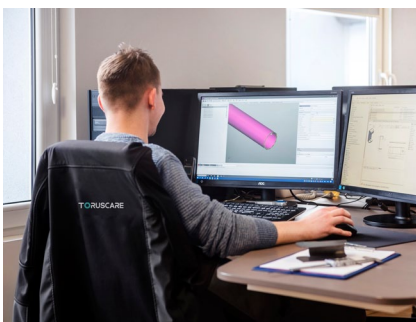

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400