



Patient Transportation Trolley

PLANO 1

We Provide Customized Overall Solutions
for Medical Equipment



PRODUCT SPECIFICATIONS

External size (LxW)	2150x830 mm	Knee-rest adjustment	0-40°
Table height adjustment	560-890 mm	Reverse trendelenburg	0-18°
Back plate adjustment	0-75°	Trendelenburg	0-18°


Crank system	1 set
Back plate air spring	1 pc
Plastic side rail	1 set
8" castors with central locking system	4 pcs
USA hydraulic pump	2 pcs
Wall bumper wheel	4 pcs
Directional wheel	1 pc
Directional pedal	2 pcs
IV pole prevision	2 pcs
IV pole	1 set
Back Rest Gas Spring	1 pc
Oxygen Cylinder Holder	1 pc
Central Locking Break Pedal	4 pcs
Transfer Mattress	1 pc
Safe Lock for Side Rail	2 pcs
Emergency Push Handle	1 pc



GUARDRAIL CAN BE FIXED AT HORIZONTAL LEVEL

Side-rail can be Fixed at Horizontal Level Fixable at horizontal level, where patient's arm can be placed on for the administration of IV or other treatment.

Proper procedure needs to be exercised to safely transfer the patient from the stretcher to the bed. An innovative side-rail design allows medical personnel to carry out this action efficiently and with greater ease. Side-rail can be laid out on a bed positioned next to the stretcher to act as a transfer board, creating zero gap for easy and quick transfer.

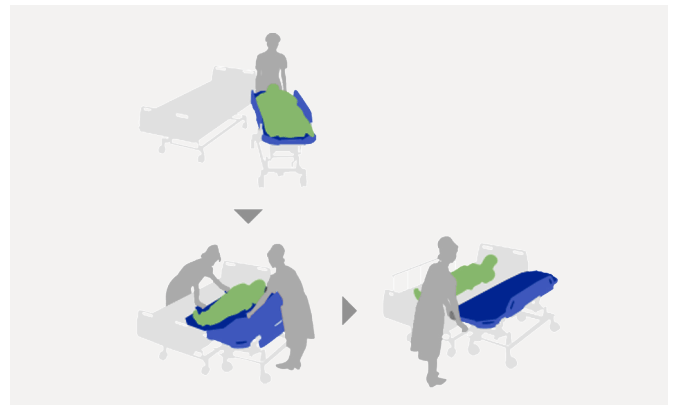

DEGREE INDICATOR

Scale on guardrail, shows the backrest degree.



T-SHAPED GAP

T-shaped design allows tubes to pass through side-rail for improved fluid drainage.



HYDRAULIC SYSTEM

American brand hydraulic pump, providing strong support for the functions of height adjustment, trendelenburg/reversed trendelenburg.

THE FIFTH WHEEL

Can be engaged to facilitate direction control, enabling users to navigate the stretcher easily in hectic and dynamic environments.



8" controlled castors.



FOLDABLE HANDLE

Equipped with foldable handle.



BUMPER

Equipped with bumper at four corners.



MATTRESS

Water proof mattress.



I.V POLE

Full SS IV pole, bearing 15kg.



POSTURE FUNCTION

Back-Rest



Knee-Rest



Bed Height Min.



Bed Height Max.




Reverse Trendelenburg




Trendelenburg



Comprehensive Hospital Design Solutions

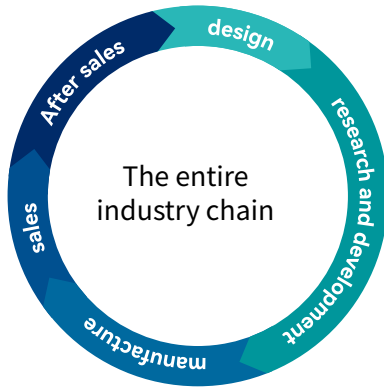


Toruscare has a complete industry chain targeting the global market.

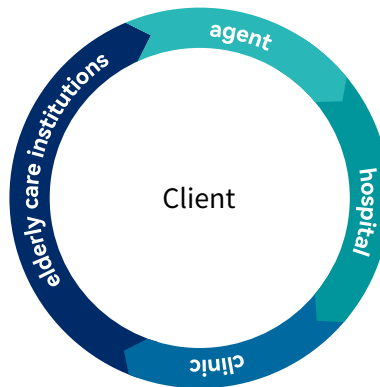

Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

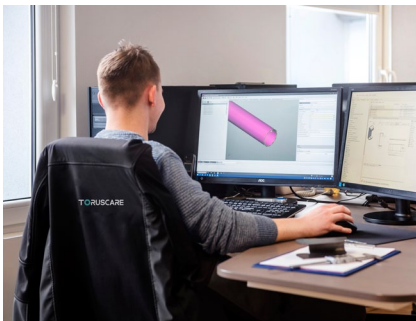

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400