



Manual Bed

GRATO 2

We Provide Customized Overall Solutions
for Medical Equipment

 www.toruscare.com

TORUSCARE

Care for All



PRODUCT SPECIFICATIONS

External Size (LxW)	2160x980mm ±10mm
Mattress Platform	1930x830 mm
Safe Working Load	200 kg
Height Adjustment	505 mm
Back-Rest Adjustment	0-80° (±10°)
Knee-Rest Adjustment	0-35° (±10°)

Manual Cranks	2pcs
5" Covered Castors	4pcs
Foldable Side Rail	1set
PP Bed Ends with Safe Lock	1set
IV Pole	1set
IV Pole Provision	4pcs
Drainage Hook	2pcs

Manual Cranks	2pcs
5" Covered Castors	4pcs
Foldable Side Rail	1set
PP Bed Ends with Safe Lock	1set
IV Pole	1set
IV Pole Provision	4pcs
Drainage Hook	2pcs


MODERN AND ELEGANT BED DESIGN

The design is modern and practical, featuring beautiful and elegant head and foot boards. The hollow protective wrap angle and barrel-type bumper wheel make the product instantly recognizable.




POSTURE FUNCTION

Back-Rest

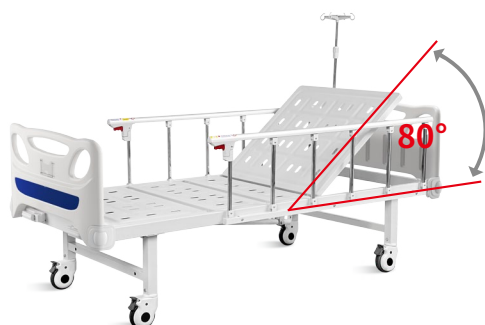

Knee-Rest

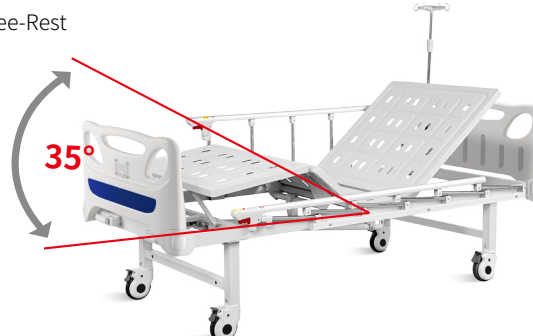

QUALITY CONTROL

Design

The design is both modern and practical. The head and foot boards feature a beautiful and elegant shape. The use of a hollow protective wrap angle, combined with a barrel-type bumper, makes the product instantly recognizable.



Steels | welding

Baosteel from Fortune Global 500.

Panasonic Robotic ensure 360° full smooth welding.



Plastic | Painting

Environmental plastic can be sterilized at a temperature of 100°C and has a tensile strength of up to 30MPa.

11 Process epoxy painting with a thickness of 0.12mm and a brightness of 60°, which is capable of resisting impacts of up to 50kg.



Accessories



Head & foot board

The head and foot boards feature "integrated protective corners" with hollow design, which are combined with barrel-type bumper wheels that have a diameter of 60mm.

This design not only minimizes the overall space occupied by the bed but also provides protection.

Additionally, the head and foot boards are effortless to disassemble and clean.

Castors

5" covered castors.

TPR material, walks 30KM with no wear and tear. Precision-designed anti-winding shell with no screws.

Pass Dynamic test: withstand a weight of 120kg, walk 30KM, and pass through 500 obstacles without any deformation or damage.



Platform

Stamping forming bed platform features a robot-welded surface that has been polished to perfection, leaving no welding scars. Bearing 200kg.



Siderail

The siderails are made of aluminium alloy horizontal bar, 6 vertical stainless steel poles and a plastic handle.



Hand Cranks

Two retractable hand cranks to control backrest and kneerest. The hand cranks remain intact after 100,000 uses.



Drainage Hooks

Injection molded drainage hooks on both sides of the bed.



IV Pole

The infusion pole is made of stainless steel and has four plastic hooks that can bear up to 15 kg of weight.



Comprehensive Hospital Design Solutions




Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

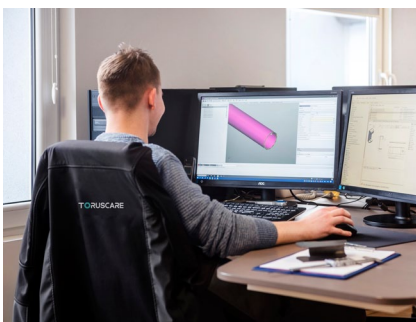

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400