



Nursing Trolley

IMPETU CL1

We Provide Customized Overall Solutions
for Medical Equipment

PRODUCT SPECIFICATIONS

Overall Dimensions	965×610×1100mm
Work Board	340×250mm
Upper Left Drawer Size	440×75mm
Lower Left Drawer Size	440×160mm
Upper Right Drawer Size	440×160mm
Lower Right Drawer Size	250×250mm

Packing Size	1300×950×650mm
Product Weight	110 kg

Drawer Divider System
Silent Directional Casters
Sealed Handle

1pc
1pc
1pc



QUALITY CONTROL

Design

The overall design and production are exquisite, the styling and color matching are reasonable, installation is simple, and the structure is stable and durable.

Material

The ABS material is an innovative engineering plastic. Antibacterial, moisture-proof, anti-corrosion, wear-resistant, non-toxic and tasteless, never rust. Washable, easy cleaning.



Drawer

Several-layer drawer design with small compartments in the middle, can accommodate multiple types of items, with powerful classification and easy and convenient access to items.



Drawer slide

The drawer uses high-strength three-section silent damping slides. When the drawer is closed to the final distance, hydraulic pressure is used to slow down the speed and reduce the impact force, creating a comfortable closing effect. Even when force is used to push the drawer, it will still gently close, ensuring perfect and quiet movement.



Drawer compartment

The interior of the drawer is equipped with transparent acrylic classified medicine grids, which can be used for classifying and placing first aid medicines. The medicine grids can be freely combined.



Working table

The trolley is equipped with a collapsible working table, which can be used for writing;

Work board size: 340*250mm;



Sealed handle

The trolley is equipped with a sealed handle, which can lock all drawers at the same time.



Castor

At the bottom, it uses 5-inch silent castors, anti-winding, high-end, silent castors, allowing for flexible movement and a high load-bearing capacity.



Comprehensive Hospital Design Solutions




Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

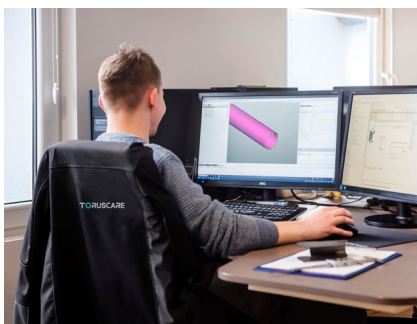

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400