

Children Bed



Children Bed

NINO S2

We Provide Customized Overall Solutions
for Medical Equipment

 www.toruscare.com

TORUSCARE

Care for All



PRODUCT SPECIFICATIONS

External size (LxW)	2030x1045 mm
Bed height	500-770±10 mm
Mattress platform	1705x900±10 mm
Back-rest adjustment	0-75° (±10°)
Knee-rest adjustment	0-35° (±10°)

Manual Cranks	3pcs
5" Hardshell Castors	4pcs
Aluminium Alloy Full Length Side Rail	1set
ABS Grace Bed Ends with Safe Lock	1set
IV Pole	1pc
IV Pole Prevision	4pcs
Drainage Hook	2pcs



Accessories

Design

Head foot board, with curve surrounded design, provide conforing feeling for patients.



Plastic

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.



Steels | Painting

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



Welding | Chassis

Panasonic Robot ensure 360° full smooth welding.



Hand Control

The Bed With Three Functions, Controlled by manual cranks.



Platform

Punching platform, smooth surface no welding scar, bearing 250KG load.



Guardrail

Sliding full side rails with distinctive design, elegant appearance; The lock bar is in conspicuous red, easy operation, sliding upper and down smoothly; Anti-pinch design, effectively protect the safety of children in use.



Head & Foot Board

Head and foot board can be reversible, extra height design. ABS Grace Bed Ends with Safe Lock, automatic lock without any manual operation on lock, on/off removable mark on hook lock.



Castors

TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.



Drainage Hooks

Mobile plastic drainage hooks.



IV pole

SS IV pole with plastic hooks, bearing 5kg.



Comprehensive Hospital Design Solutions




Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

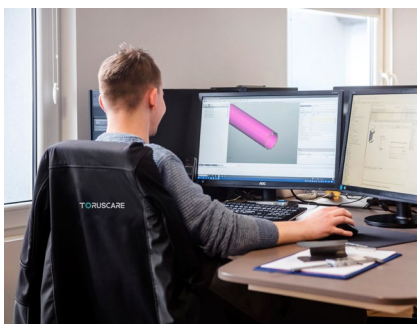

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400