



Electric Dialysis Chair

Electric Dialysis Chair

QUILAR 2

We Provide Customized Overall Solutions
for Medical Equipment

PRODUCT SPECIFICATIONS

Size	1910×580mm
Back	780×580mm
Seat	480×580mm
Leg	530×580mm
Seat height	550-700mm
Armrest	500×180mm

With CPR
Hand controller
Power cord
3" silent casterPP Bed Ends with Safe Lock
Small dining table


Mattress thickness	70mm
Back rest lifting	-8°-70°
Leg rest lifting	-50°-0°
Armrest swing out	0°-270°
Foot board adjustable distance	150mm
Voltage and frequency	100-240 V · 50/60 Hz
Plug as required (hospital grade)	

1pc
1pc
4set
1pc

QUALITY CONTROL

Material

Baosteel from Fortune Global 500.

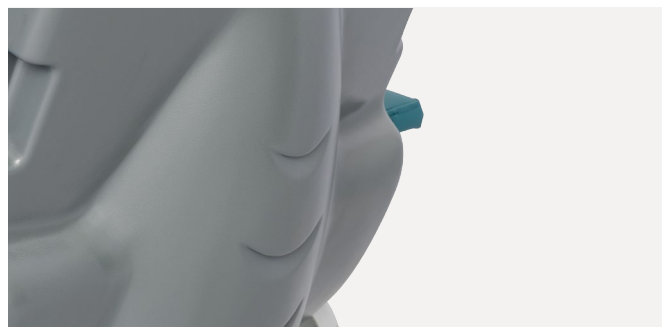

Plastic

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



Welding

Panasonic Robotic ensure 360° full smooth welding.



Electric Motor

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.



Chair Cover

Optional genuine leather or supebric cover, with artisanal making, comfortable hand feeling.



Table Board

Adjustable ABS tray table, with antiskid convex edge, tank for cup and cellphone. USB/5V optional. Very convenient for donor.



Handrail

Upgraded armrest, can be right-left, front-back adjustable. Donor can control universal position.

Armrest is R120° concave curve, avoid donor's arm sliding,ensure the donor's safety and comfortable.



Foot Board

150mm foot board adjustable distance.



Ultrafiber leather color table



* Note: the above is the physical scan, there may be a certain color difference in kind.

Comprehensive Hospital Design Solutions




Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

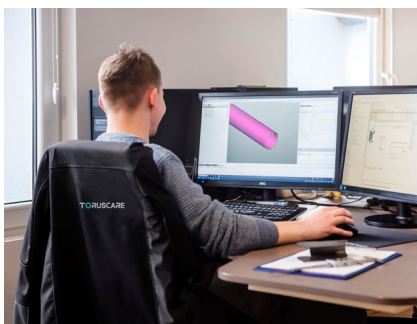

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400