



Electric Homecare Bed

SABER 5

We Provide Customized Overall Solutions
for Medical Equipment



PRODUCT SPECIFICATIONS

External size (LxW)	2060×1070mm	Trendelenburg	0-15°
Mattress platform	2000×900mm	Lowest & Highest Position	300-800mm
Back-rest Adjustment	0-68°	Safe working load	215kg
Leg-rest Adjustment	0-28°	Mattress size	L2000×900mm/800mm
Reverse Trendelenburg	0-17°		

Linear Motor	4pcs
Wooden side rails	1set
Head & Foot Board	1set
Castors	4pcs
Remote controller	1pc
Emergency power system	Optional
Trapeze Pole	Optional



KEY FEATURES

- Guardrails divided in a 40/60 ratio
- Modern pantograph design
- Various wood colors to choose from
- Wide range of height adjustment (lowest position of the bed – 30cm)
- Head and foot boards independent of the bed structure
- Maximum working load: 215kg
- 2-year warranty

SAFETY LABELS

multiple warning labels on the bed body to prevent pinching and potential risks and injuries.



OVERLOAD PROTECTION

The motors are equipped with electric protection against overload, therefore the bed will automatically stop when a motor is overloaded. Once the overload is removed, the bed can continue to work without interference.



MATTRESS PLATFORM

The bed board is divided into 4 parts, the back and legs can be folded up separately, and the size division is ergonomic.



HEAD & FOOT BOARD

The head and foot boards are made of laminate and solid beech wood, providing strong stability while being lightweight.



SUPER LOW POSITION

Lowest height 300mm from the ground, to prevent accidental falls from the bed.

The bed can be adjusted to the right height, which can reduce efforts to get in and out of the bed, wash feet, give daily nursing care and other movements.



SIDERAIL

4-piece guardrail, three-profile structure, one-touch folding, curved shape and covered seam design, increase space and improve comfort.




REMOTE CONTROLLER

The 10-button handset is equipped with a locking function to prevent accidental touches.



ADJUST DESIGN TO INDIVIDUAL NEEDS

Wide range of wood colours

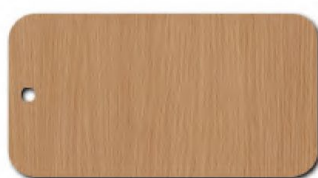

Beech



Light oak




POSTURE FUNCTION

Back-Rest



Knee-Rest

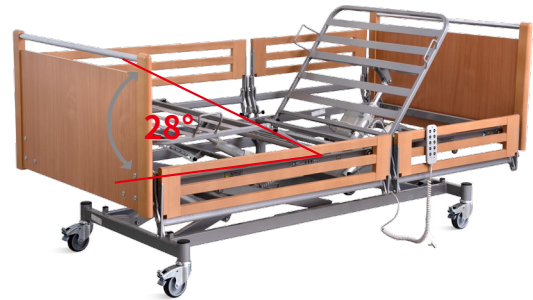

Bed Height Min.



Bed Height Max.




Reverse Trendelenburg



Trendelenburg



Comprehensive Hospital Design Solutions

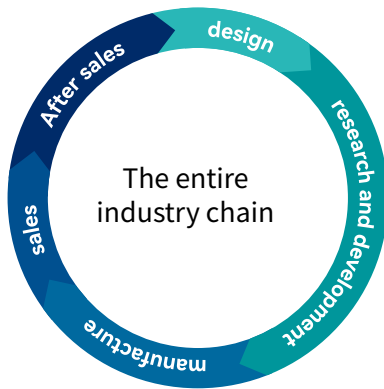


Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

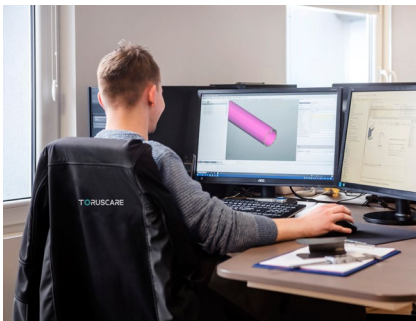

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400