



Electric Operating Table

CURAR 4

We Provide Customized Overall Solutions
for Medical Equipment

PRODUCT SPECIFICATIONS

External Size (LxWxH)	2150x555x(700-1000) mm
Reverse Trendelenburg	0-25°
Trendelenburg	0-25°
Lateral tilt left/right	0-25°
Head plate up/down	45°/90°
Head plate Height adjustment	70 mm
Back plate up/down	70°/22°
Leg plate up/down	15°/90°

Leg plate outwards	90°
Kidney bridge elevation	0-110 mm
Horizontal sliding	0-350 mm
Safe Working Load	275 kg
X-rays translucent	YES
C-amm compatible	YES
Multi-functional orthopedic frame	YES

Bed frame
Mattress
Anesthesia screen frame
Shoulder holder
Arm holder
Leg support
Foot support
Hand controller
Sliding block
Spanner

1pc
1set
1pc
2pcs
2pcs
2pcs
2pcs
1pc
8pcs
5pcs

Material

The base material for the table is high thickness 304 stainless steel. seamless welding by robot, with modern appearance and more stable.




Mattress

Polyurethane foam mattress, waterproof and seamless, avoid liquid in. The connection between table and mattress is buckle joint, easy to disassemble and install, easy to clean.



Head plate

Head plate is with imported gas spring for easy control.

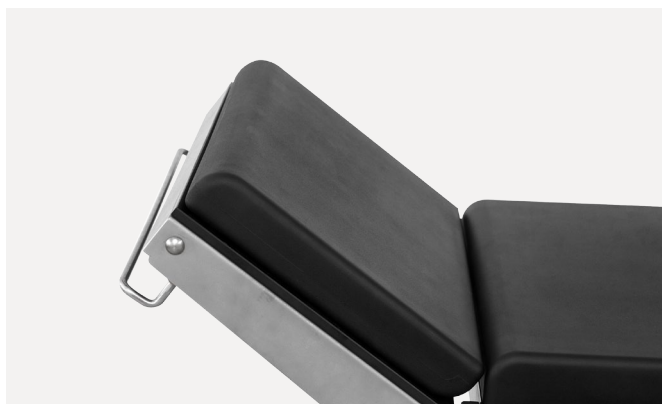

Kidney bridge

Manual kidney bridge, operate by head side rocker, more convenience.



Horizontal sliding

Extra horizontal sliding, distance 350mm high strengthen phenolic resin table surface, C-arm and X-ray compatible.



Dual control system

All functions controlled by both hand controller or control panel on table base column. With action LCD display on hand controller, easy to operate.

With control panel on table base column, in case of emergency hand controller failure.



Leg plate

The up and down folding actions of the leg plate are all completed by electric.

The leg plate is removable, effortless and smooth.



Castors

HAION brand castors, each load bearing is 250kg, smooth and noiseless movement, with longer service life.

Pedal brake, easy to move and lock. Unique switch design.



Battery


Battery could support about 80-120 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table runing in the event of power failure.



Armrest



Anesthesia Frame

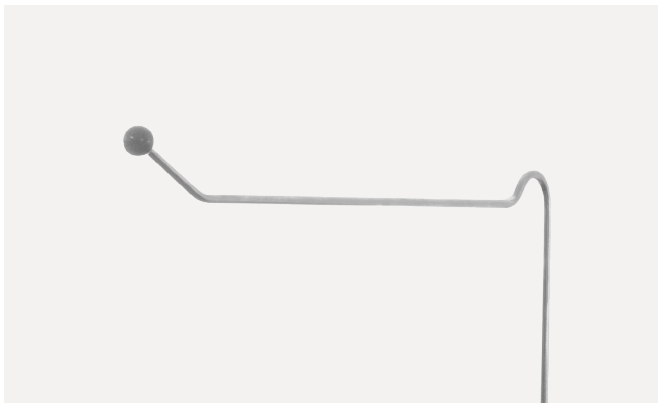

Leg Crutches



Shoulder support



Body supporter



POSTURE FUNCTION

Lowest position of table top

700 mm




Highest position of table top

1000 mm




Back plate up

70°

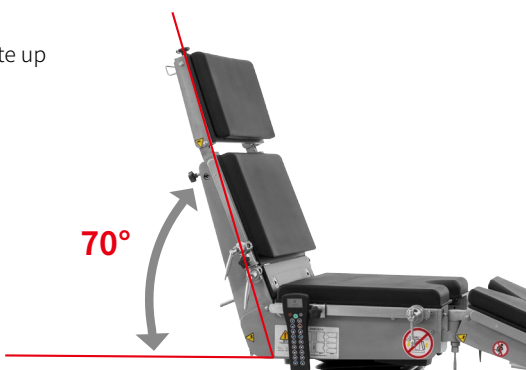

Back plate down

22°

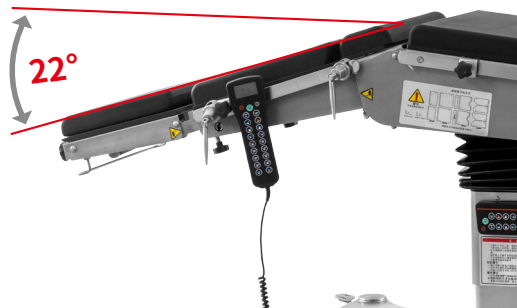


Reverse Trendelenburg

25°



Trendelenburg

25°

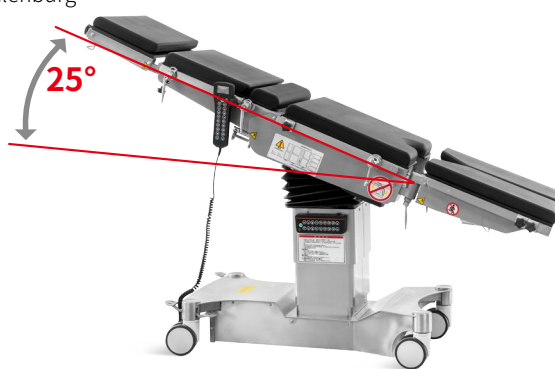

Lateral Tilt Right

25°



Lateral Tilt Left

25°



Head plate upward



Head plate downward



Leg section evagination

0-90°



Leg plate down

0-90°



Head plate Height adjustment

0-70mm



Kidney bridge elevation

0-110mm



Horizontal sliding

0-350mm




Comprehensive Hospital Design Solutions




Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.

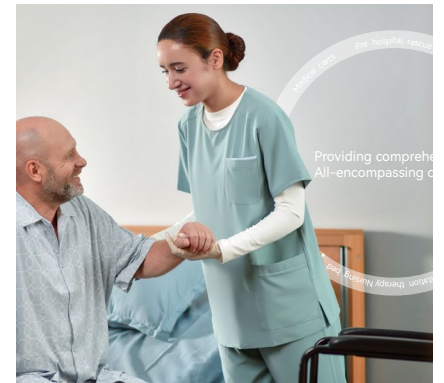

Toruscare provides services for full coverage of target audiences.

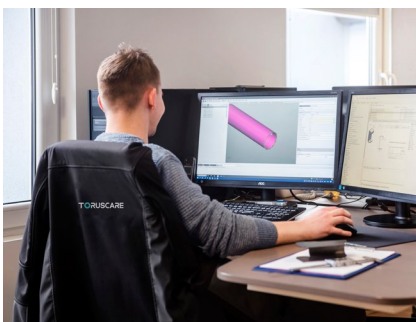

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400