

Electric Operating Table



Electric Operating Table

CURAR 7

We Provide Customized Overall Solutions
for Medical Equipment

PRODUCT SPECIFICATIONS

External Size (LxWxH)	2080x550x(700-1000) mm
Reverse Trendelenburg	0-25°
Trendelenburg	0-25°
Lateral tilt left/right	0°-25°/0°-25°
Head plate up/down	0°-50°/0°-90°
Back plate up/down	0°-80°/0°-20°
Leg plate up/down/outwards	0°-20°/0°-90°/0°-90°

Horizontal sliding	0-300 mm
Kidney bridge elevation	0-140 mm
Safe Working Load	250 kg

Shoulder support
Arm support
Lumbar support
Leg support
Foot support
Sliding block
Spanner
Hand controller
Mattress
Anesthesia screen frame

2pcs
2pcs
2pcs
2pcs
2pcs
8pcs
5pcs
1pc
1set
1pc

Head Plate

Uses a gas springsupport, ensuring quiet operation.
Can stop at any angle during operation.



Leg Plate

Utilizes a vertebra link with gas spring support, allowing for one-handed and effortless operation.
Can stop at any angle; Outward rotation requires only a slight loosening of the locking axis to achieve leg board abduction.



Motor Cable

Concealed motor cable, keeping the underside of the tabletop clean and unobstructed for clear imaging.




Hand Controller

Connection cable is located under the base, keeping the underside of the tabletop clean and unobstructed for clear imaging. Features a one-click plug-and-play design for easy replacement and maintenance. Automatically shuts off after 15 minutes of inactivity. Includes an emergency stop button and battery level display.



Mattress

High frequency heat sealing technology.
Large area of velcro coverage provides strong adhesion.

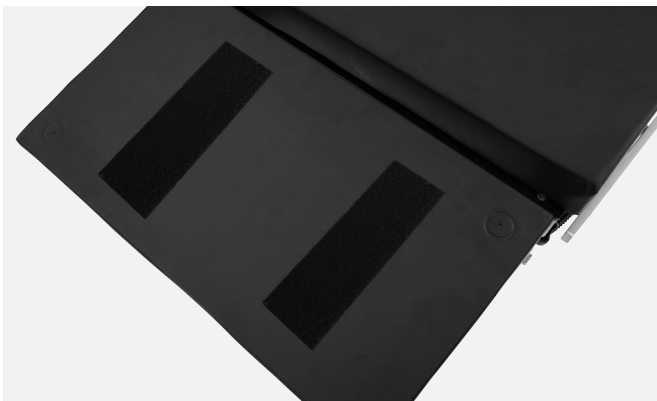


Hand Controller Storage

Can be hung at any position on the headboard without affecting clinical use.



Battery

Battery could support about 80-120 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table running in the event of power failure.

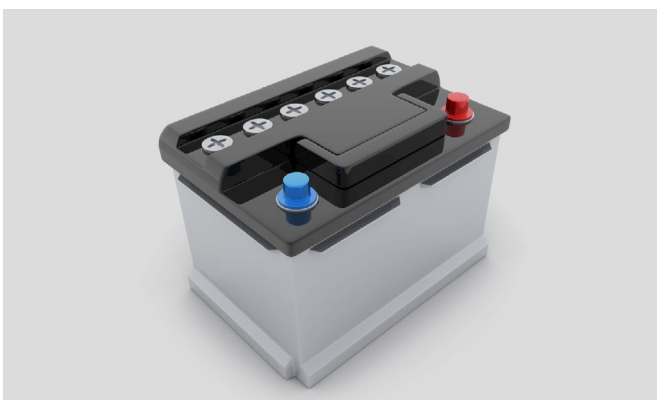

Kidney Bridge

Manual kidney bridge, operate by head side rocker, more convenience.



Anesthesia Screen Frame

Equipped with a stainless steel anesthesia screen frame for hanging anesthesia curtains, used to shield and cover the patient's head during use.



Shoulder Support

Polyurethane support frame, adjustable in angle, used to stabilize the patient's body.



Body Support

Polyurethane support frame, adjustable in angle, used to stabilize the patient's body.




Arm Rest

Polyurethane widened armrest board with attached straps providing stable performance and versatile applications.



Leg Rest

Manual multi-angle adjustment for easy use.

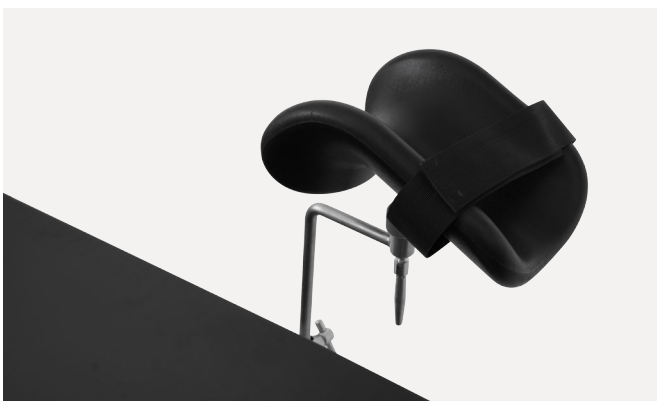

Foot Plate

Be used to support patient's feet.



POSTURE FUNCTION

Lowest position of table top

700 mm



Highest position of table top

1000 mm



Back plate up

0°-80°




Back plate down

0°-20°




Reverse Trendelenburg

25°

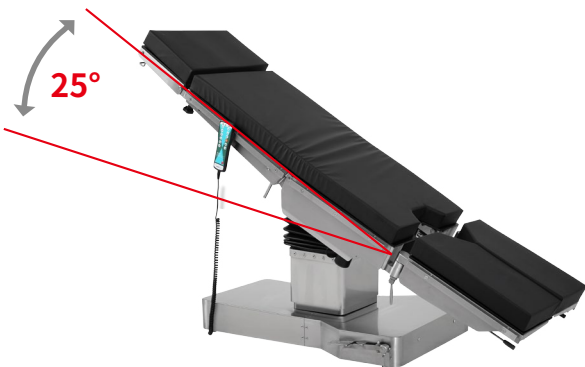

Trendelenburg

0-25°

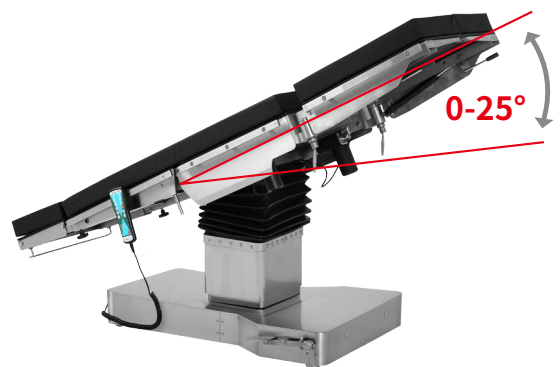

Left Lateral Tilt

25°




Right Lateral Tilt


25°



Leg plate up

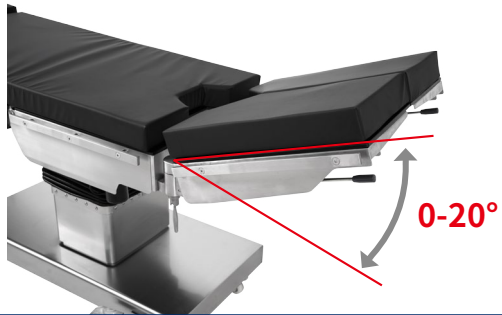

Head plate downward

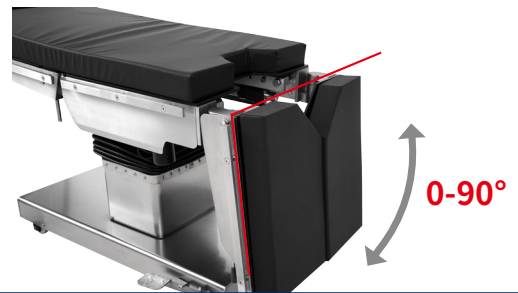

Leg plate outwards

0-90°




Kidney bridge elevation

0-140mm



Head plate up

0-50°

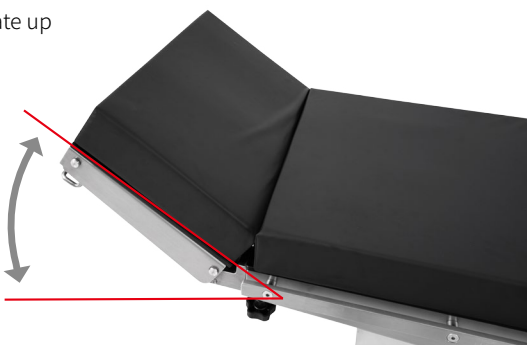

Head plate down

0-90°



Horizontal sliding

0-300mm



Comprehensive Hospital Design Solutions




Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

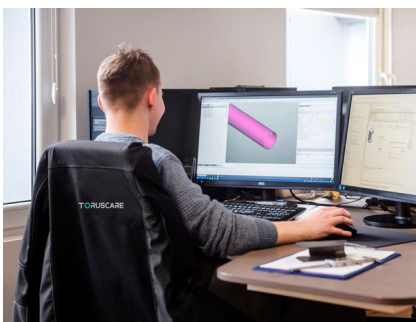

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400