



Rehabilitation Therapy Nursing Bed Medical Carts Pre Hospital Rescue Supporting Facilities Operating Room Gynaecology

Electric Children Bed


NINO S

We Provide Customized Overall Solutions
for Medical Equipment

 www.toruscare.com

TORUSCARE

Care for All

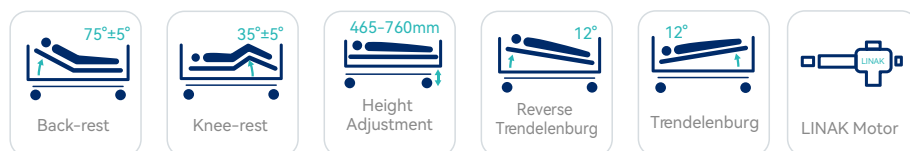


PRODUCT SPECIFICATIONS


External Size (LxW)	1975x880 mm	Reverse Trendelenburg	0-12°
Mattress Platform	1730x780 mm	Trendelenburg	0-12°
Safe Working Load	200 kg		
Height Adjustment	465-760 mm		
Back-Rest Adjustment	0-75° (±5°)		
Knee-Rest Adjustment	0-35° (±5°)		
Linak Motor		4pcs	
Stainless steel side rail with extra protection		1set	
PP Bed Ends with Safe Lock		1set	
PP Bed Platform, removeable, washable		5pcs	
5" doubleside castors with central locking system		4pcs	
IV Pole		1set	
IV Pole Provision		4pcs	
Drainage Hook		2pcs	
Trapeze pole		Optional	
Directional wheels function		Optional	

SIDE RAIL MADE FROM STAINLESS STEEL AND TAKE 12-STAFF DESIGN

Our beds feature a unique aseptic design tailored to meet hospital infection control requirements. The seamless integration of head, tail, and bed surfaces supports swift and thorough disinfection processes.

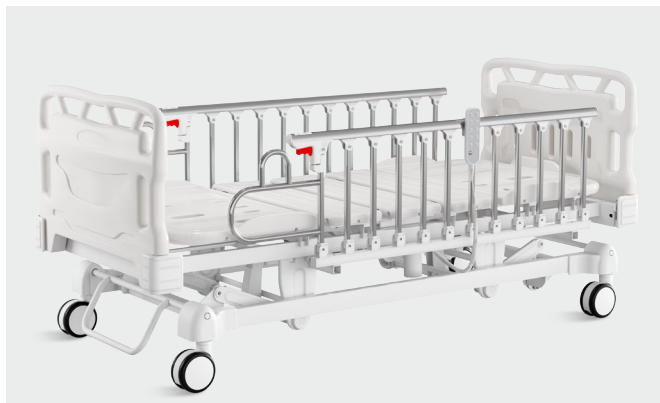

The stainless steel side rails are designed with a 12-staff configuration, providing robust support. They include stainless steel protection to facilitate easy access for patients getting on and off the bed.

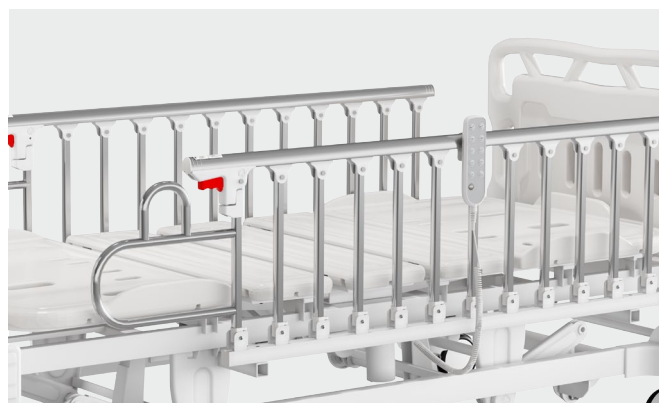

QUALITY CONTROL

Design

This Bed uses unique design to meet hospital strict gremfree request. Head & Foot board and platform are all seamless, which supports effective sterilization.



Steels | welding

Baosteel from Fortune Global 500.

Panasonic Robotic ensure 360° full smooth welding.



Plastic | Painting

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

11 Process epoxy painting, ASTM testing anti-baterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



Electric Motor

American motor, with UL/ROHS/EN standard. Noise under 50dB, with lifetime 20,000 times movement.



Accessories

Platform

Five PP material bed platforms are washable and easily movable.



Castors | Brake system

Castors with TPR tire, no worn-out after running 30KM, united forming without bolts.

Pass Dynamic test: bearing 120kg run 30KM and pass obstacles 500time.

Brake uses 301 stainless steel. Casters' brakes keep qualified after 10,000 times test.



Head & Foot Board


Modern design PP head&foot board.

Easy to be cleaned and installed.




Handle Controller

Handle controller contains flexible spiral cable for easy mobility and operation in optimal position. Patient and nurse can control all necessary functions.

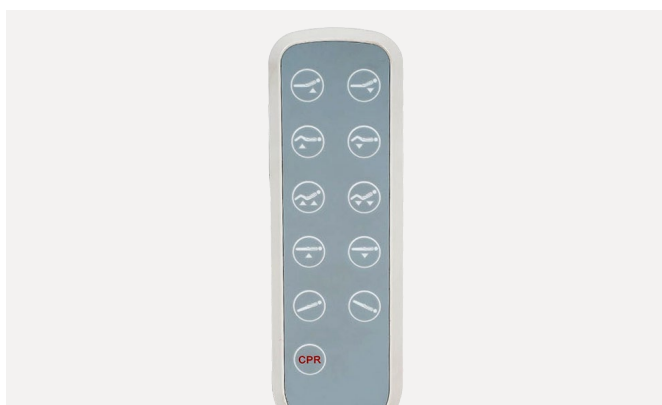

Siderail

Side rail made from stainless steel and take 12-staff design. And it with stainless steel protection for patients can get on or off bed easily.

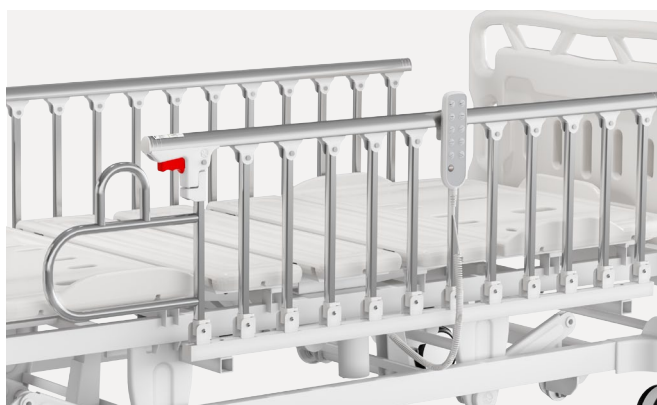

Drainage Hooks

Mobile plastic drainage hooks.



IV Pole

Full SS IV pole, bearing 2.5kg per hook.




POSTURE FUNCTION


Back-Rest




Knee-Rest

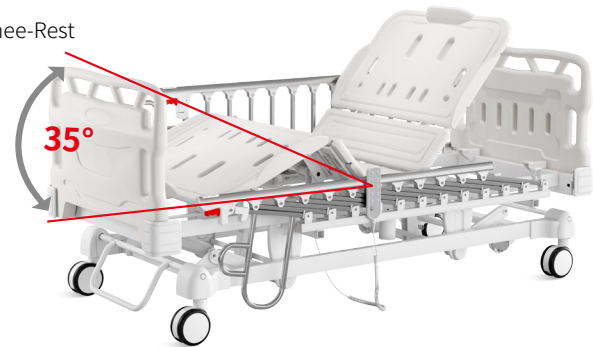


Bed Height Min.

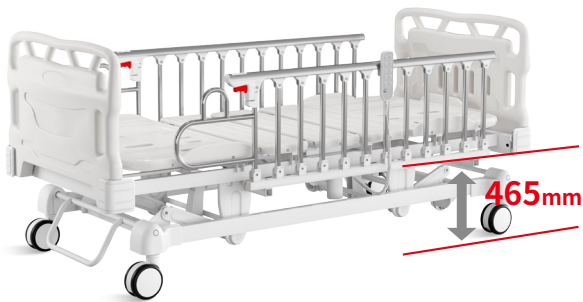


Bed Height Max.

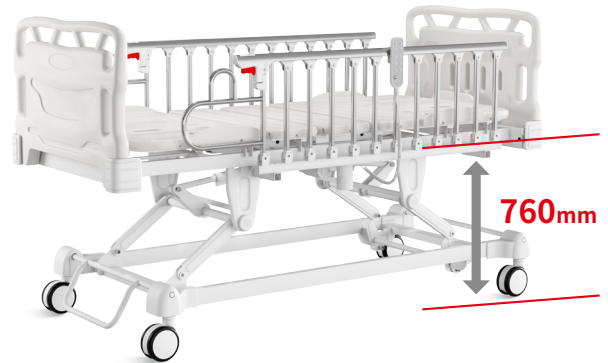

Reverse Trendelenburg

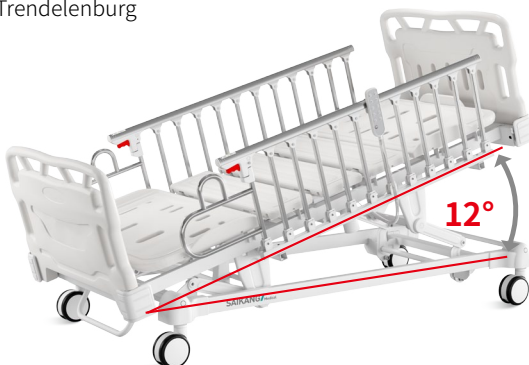

Trendelenburg

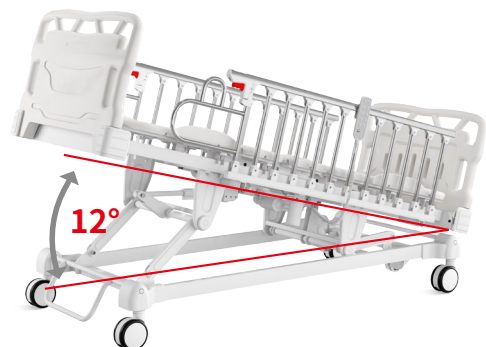

Comprehensive Hospital Design Solutions




Toruscare has a complete industry chain targeting the global market.



Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

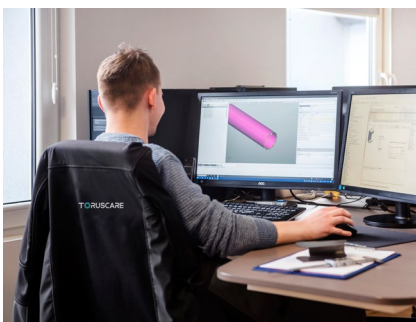

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400