

Care for All



Patient Transportation Trolley


PLANO 2

We Provide Customized Overall Solutions
for Medical Equipment

 www.toruscare.com

TORUSCARE

Care for All

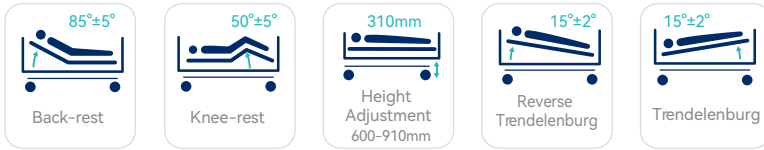

PRODUCT SPECIFICATIONS

External size (LxW)	2200x770 mm	Reverse trendelenburg	0-15°±2°
Table height adjustment	600-910 mm	Trendelenburg	0-15°±2°
Back plate adjustment	0-85°±5°	Static load	250kg
Knee-rest adjustment	0-50°±5°	Dynamic load	135kg

80mm PU Mattress	1pc
8" Castors, Central Locking	4pcs
Hydraulic Pumps	2pcs
IV Pole	1set
Back Angle Monitor	2pcs
Integrated Wall Bumper Wheels	4pcs
Whole Body X-ray with Cassette Holder	1set
Patient Restraints	2pcs
PP Bottom Cover	1pc
5-speed Foldable Guardrail	1set
Foot Pedals	1set
5th Directional Wheel	1pc

OPTIONAL

Oxygen Tanks Holder
Paper Holder
Foot-Board & Monitor Holder
Monitor Holder
Debris Basket
Weighting System
Pluggable Mobile IV pole



SUITABLE FOR PATIENT TRANSPORT IN HEALTHCARE FACILITIES

Providing comfort for both patients and caregivers.

Composed of the trolley body, two hydraulic cylinders, the back gas spring mechanism, the leg rocker mechanism and the foot pedal, these systems which can realize adjustment of trolley's position.



PLASTIC | PAINTING

Environmentally friendly plastic is non-toxic and tasteless, and can be sterilized at temperatures above 100°C.

Vivid color.



HYDRAULIC SYSTEM

Two hydraulic cylinders provide strong support for the lifting and trendelenburg/reversed trendelenburg of the trolley.




GUARDRAIL


Foldable high-strength 5-gears guardrail to maintain patient safety;
Guardrail use method: Open the guardrail switch, push guardrail to the direction of the leg; Pull up the guardrail, lock the guardrail with the guardrail switch, and keep the guardrail standing.

THE FIFTH WHEEL

The diameter of the fifth wheel is 125mm, which can control the trolley go straight when it is lowered, and the trolley can go in all directions when it is retracted.

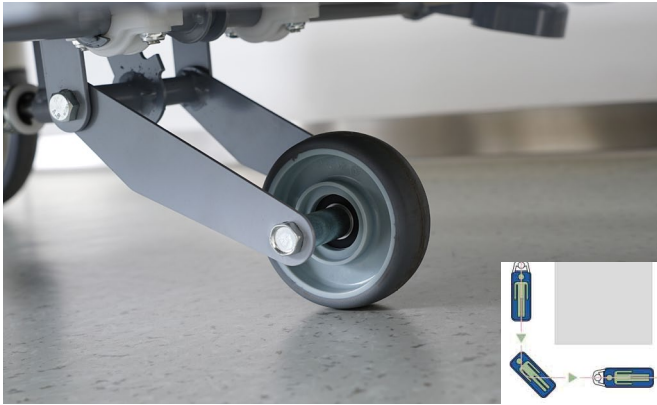

Four 200mm diameter central locking castors, the dust-proof can prevent hair into the caster to cause the inconvenience of movement, both sides with central locking casters, one pedal control all 4 castors.

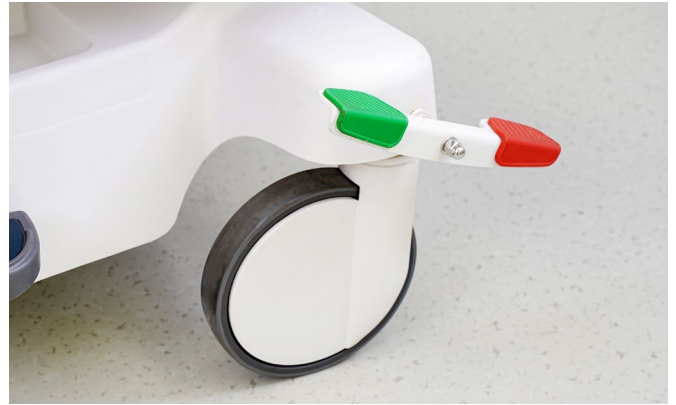

PLATFORM

Phenolic resin bed board with cassette holder meets radiographic requirements.



X-RAY HOLDER

X-ray holder on the back.



BASE COVER

PP base cover which can place items and oxygen cylinder.



BUMPER

4 built-in removable and durable bumper wheels, diameter 75mm±1mm.



Mattress stopper

Both the head and the end of the bed are equipped with mattress stopper to prevent the mattress from sliding.




Hand bars

Head of equipped with 2 hand bars for convenience of medical staff.



Body strip

2 body straps.

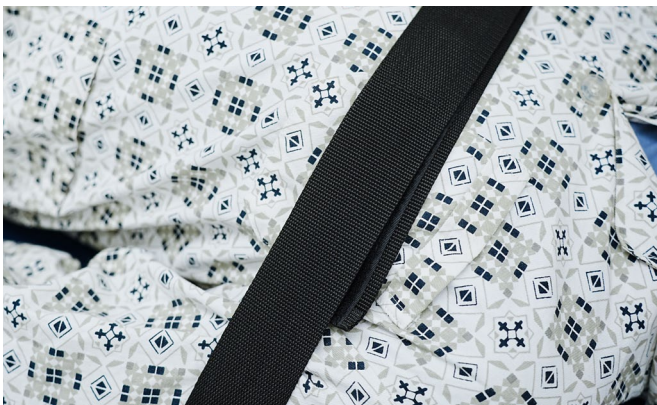

Angle indicator

2 angle indicators inserted near the head and back of the trolley.



Mattress

8cm PU mattress.



IV Pole

IV Pole provision in the head part, equipped with retractable stainless steel with 4 hooks.



OPTIONS ACCESSORIES

OXYGEN TANKS HOLDER

Vertical oxygen tank holder, which can be placed in 4 corners(removable).



PAPER HOLDER

Head roll holder (removable).



FOOT-BOARD&MONITOR HOLDER

Removable insert plate & plastic instrument platform.




MONITOR HOLDER

Phenolic resin instrument platform.



DEBRIS BASKET

Hanging basket(end part).

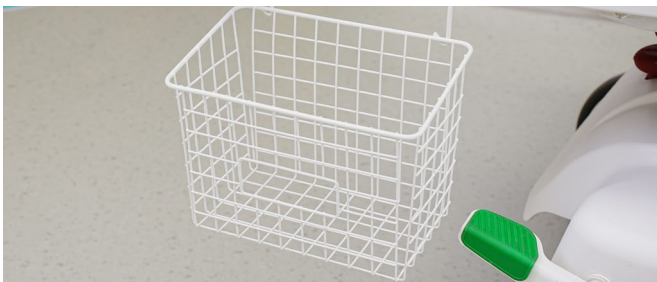

WEIGHTING SYSTEM

Weighing system (LCD display).



IV POLE

Removable mobile IV Pole.



POSTURE FUNCTION

Back-Rest



Knee-Rest



Bed Height Min.



Bed Height Max.




Reverse Trendelenburg




Trendelenburg



Comprehensive Hospital Design Solutions

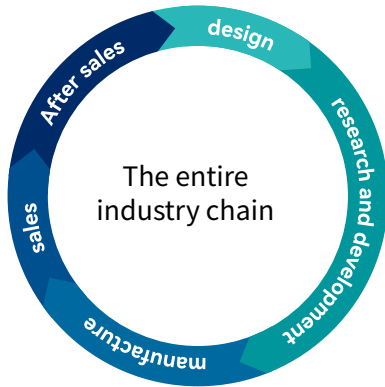


Toruscare has a complete industry chain targeting the global market.

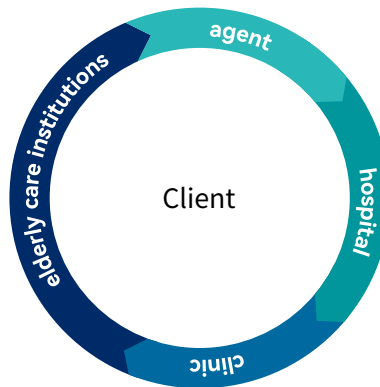

Toruscare provides overall healthcare solutions.



Toruscare provides services for full coverage of target audiences.

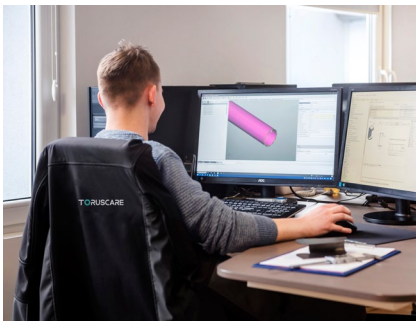

Design and Development Excellence

At Toruscare, our strength lies in our robust design and development capabilities. Our dedicated team of engineers and designers continually push the boundaries of innovation, ensuring that our products are at the forefront of medical technology. From conceptualization to prototyping, we meticulously craft each product to meet the highest standards of functionality, durability, and ease of use.



Uncompromising Quality Control

Toruscare is deeply committed to quality, implementing stringent control measures at every manufacturing stage. Our state-of-the-art facilities, equipped with advanced testing and inspection tools, ensure each product meets rigorous standards. This dedication to excellence has earned the trust and loyalty of healthcare providers worldwide.



Global Reach, Local Touch

As a global company, Toruscare understands the diverse needs of different markets. We pride ourselves on our ability to deliver customized solutions that cater to the unique requirements of healthcare providers across the world. Our extensive distribution network and strategic partnerships ensure timely delivery and exceptional after-sales support, making us a trusted partner in the global healthcare community.



Care for All

TORUSCARE

Toruscare Sp. z o.o

Adress: UL. SPACEROWA 1C, 41-253 CZELADŹ, POLAND

Email: export@toruscare.com

Phone: +351 910200400